LED ROVER

Instructions for Use

ST Technologies
6018 Bowdendale Avenue
Jacksonville, FL 32216 USA
Customer Service: 904 208 2291
Toll Free 877 814 2237

RMS UK, Ltd.
28 Trinity Road
Nailsea, Somerset BS48 4NU
United Kingdom
TEL: 01275 858891

TABLE OF CONTENTS

1. INTRODUCTION
2. INTENDED USE/INDICATIONS FOR USE
3. WARNINGS/LABELING
4. SPECIFICATIONS
5. OPERATION
6. BATTERY INFORMATION
7. CLEANING/STERILIZATION
8. WARRANTY
9. CHART OF MEDICAL DEVICE STMBOLS

1. INTRODUCTION

The LED Rover is a Light Emitting Diode (LED) and battery unit to couple onto the light post of commercially available endoscopes having an ACMI or Wolf light post. The LED Rover is made to take the place of conventional light cable and remote light source. It couples to the endoscope in the same fashion as a conventional light cable by press fitting onto the endoscope's light post. The LED Rover is powered by two lithium 3V, CR123 batteries (or interchangeable equivalents).

The batteries power an LED via a current converter that supplies a constant current to the LED. The light is directed into the light post of the scope. The light is then carried by fiber optics internal to the scope to illuminate the target.

2. INTENDED USE / INDICATIONS FOR USE

The LED Rover is used to illuminate the site of surgery during minimally invasive surgical procedures in arthroscopy (orthopedic surgery), laparoscopy (general and gynecological surgery) and in Endoscopy (general, gastroenterological and ENT surgery). The light is transmitted from the light source through a scope.

3. WARNINGS / LABELING

3.1 WARNINGS

Caution Battery powered equipment is considered unsuitable for use in the presence of flammable anesthetics.

Caution Federal law restricts this device to sale by or on the order of a licensed practitioner.

Caution The LED Rover is provided **non-sterile**. If used for sterile procedures, the unit must be sterilized or draped with commercially available surgical drapes.

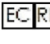
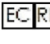
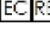
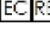
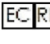
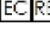
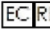
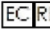
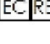
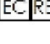
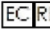
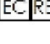
 **Caution** The LED Rover will become warm to the touch but even prolonged contact during use will not cause injury.

 **Caution** Avoid looking directly at the LED inside the port or pointing it directly at anyone's eyes. The LED is painfully bright, can momentarily blind a person and will cause lingering ghost images.

 **Caution** Handle damaged or leaking batteries with care. Should you come in contact with the electrolyte, wash the exposed area with soap and water. Should the eye come in contact with the electrolyte, flush the eye with water for 15 minutes and seek medical attention. Do not incinerate discarded batteries. Follow local EPA requirements for proper disposal.

3.2 LABELING

The LED Rover will be wrapped in bubble wrap and placed in the white box labeled as shown below.

<p align="center">CUDA[®] <i>SURGICAL</i></p> <p align="center">ROVER III LED LIGHTSOURCE, ACMI PORT</p> <p align="center">Caution: Federal law restricts this device to sale by or on the order of a licensed practitioner.</p> <table border="0"> <tr> <td></td> <td>LLS-3LBA</td> <td rowspan="2">R_x</td> </tr> <tr> <td></td> <td>WOXXXXX</td> </tr> <tr> <td></td> <td>ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.</td> <td rowspan="2">CE</td> </tr> <tr> <td></td> <td>20YY</td> </tr> <tr> <td></td> <td>RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom</td> <td></td> </tr> </table>		LLS-3LBA	R_x		WOXXXXX		ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.	CE		20YY		RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom		<p align="center">CUDA[®] <i>SURGICAL</i></p> <p align="center">ROVER III LED LIGHTSOURCE, ACMI PORT RECHARGEABLE BATTERIES</p> <p align="center">Caution: Federal law restricts this device to sale by or on the order of a licensed practitioner.</p> <table border="0"> <tr> <td></td> <td>LLS-3LBA-R</td> <td rowspan="2">R_x</td> </tr> <tr> <td></td> <td>WOXXXXX</td> </tr> <tr> <td></td> <td>ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.</td> <td rowspan="2">CE</td> </tr> <tr> <td></td> <td>20YY</td> </tr> <tr> <td></td> <td>RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom</td> <td></td> </tr> </table>		LLS-3LBA-R	R_x		WOXXXXX		ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.	CE		20YY		RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom	
	LLS-3LBA	R_x																									
	WOXXXXX																										
	ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.	CE																									
	20YY																										
	RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom																										
	LLS-3LBA-R	R_x																									
	WOXXXXX																										
	ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.	CE																									
	20YY																										
	RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom																										
<p align="center">CUDA[®] <i>SURGICAL</i></p> <p align="center">ROVER III LED LIGHTSOURCE, WOLF PORT</p> <p align="center">Caution: Federal law restricts this device to sale by or on the order of a licensed practitioner.</p> <table border="0"> <tr> <td></td> <td>LLS-3LBW</td> <td rowspan="2">R_x</td> </tr> <tr> <td></td> <td>WOXXXXX</td> </tr> <tr> <td></td> <td>ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.</td> <td rowspan="2">CE</td> </tr> <tr> <td></td> <td>20YY</td> </tr> <tr> <td></td> <td>RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom</td> <td></td> </tr> </table>		LLS-3LBW	R_x		WOXXXXX		ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.	CE		20YY		RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom		<p align="center">CUDA[®] <i>SURGICAL</i></p> <p align="center">ROVER III LED LIGHTSOURCE, WOLF PORT RECHARGEABLE BATTERIES</p> <p align="center">Caution: Federal law restricts this device to sale by or on the order of a licensed practitioner.</p> <table border="0"> <tr> <td></td> <td>LLS-3LBW-R</td> <td rowspan="2">R_x</td> </tr> <tr> <td></td> <td>WOXXXXX</td> </tr> <tr> <td></td> <td>ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.</td> <td rowspan="2">CE</td> </tr> <tr> <td></td> <td>20YY</td> </tr> <tr> <td></td> <td>RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom</td> <td></td> </tr> </table>		LLS-3LBW-R	R_x		WOXXXXX		ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.	CE		20YY		RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom	
	LLS-3LBW	R_x																									
	WOXXXXX																										
	ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.	CE																									
	20YY																										
	RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom																										
	LLS-3LBW-R	R_x																									
	WOXXXXX																										
	ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.	CE																									
	20YY																										
	RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom																										

The device has no known contraindications.

LIT102 CUDA[®] Surgical
(English)

Rev. O
Date of Revision: 07/22/14
Page 4 of 40

4. SPECIFICATIONS

Item	Specification
Lamp Type	LED
Color temperature	6500° K (nominal)
Lamp life	50,000 hrs. (nominal)
Battery	Lithium CR123, 3Volts (2X) or Lithium Ion RCR123A 3V (2X)
Current converter output	611 milliamps (nominal)
Discharge time	110 minutes (typical, actual endurance may vary)
Model	LED Rover
Dimensions	5.1 inches long X 1.25 inches diameter at LED housing
Weight	3.7 oz.
Electrical Protection	Class I, Type BF
Medical Classification	FDA Class II Europe Class I

5. OPERATION

The LED Rover is operated by a pushbutton switch inside a protective silicone rubber boot at the back end of the handle. It features a momentary ON function with partial depression of the switch, and clicks to stay ON with full depression. Depress the switch fully again to turn it off.

6. BATTERY INFORMATION

6.1 DISPOSABLE BATTERY MODEL BATTERY REPLACEMENT

Unscrew the headpiece that houses the LED. (Do not remove the switch; the batteries cannot fit through the switch end.) Remove the old batteries and insert the new ones. Observe correct polarity. The positive (+) end of both batteries must point toward the headpiece. The positive end of the battery has a protruding nipple, the negative (-) end is flat.

NOTE: *Because discharge characteristics may vary between designations or between used and unused batteries, always replace batteries in pairs with pairs of identical new batteries.*

Replace disposable 3 volt Lithium CR123 primary batteries only with similar size and voltage batteries.

Storage: Store in a dry area at 40° F – 85° F.

6.2 RECHARGEABLE BATTERY MODEL

The rechargeable model lightsource differs from the standard only by the use of rechargeable batteries. The rechargeable batteries provided are a high capacity “RCR123A” 3 volt Lithium-ion battery. The charger is compatible only with the 3 volt high capacity rechargeable Lithium-ion batteries. The 3 volt high capacity batteries provided are compatible only with chargers specifically designed for them. Do not put other types of batteries in the charger or attempt to charge the batteries with another type of charger.

A pair of fully-discharged batteries will take approximately three hours to reach a full charge. The indicator light on the charging cradle will show green when the power converter powers the empty cradle. The indicator is red/orange when charging batteries and green when batteries are fully charged. Batteries may be “topped off” to full charge even if they are not discharged. Batteries may be left in the charging cradle when charged. (If the cradle is not powered, the indicator will show green if a battery with even a very low charge is inserted; the cradle must be powered for the indicator to be true.) It is advised to keep pairs of batteries together in use and in charging so that they both have approximately the same state of charge. Damage can result from mixing batteries of different charge strengths either in use or in charging. Avoid directly short circuiting a battery. The power converter operates on 100~240 VAC, 50~60Hz @ 0.5A.

Fully charged batteries will operate the lightsource at full brightness for approximately 1½ hours, followed by approximately ½ hour of usable gradually diminishing light.

7. CLEANING / STERILIZATION

7.1 CLEANING

The LED Rover is a high quality optical device. Care for it appropriately.

A soft brush with mild detergent may be used to clean the device. The lens of the LED inside the optic port may be wiped with a swab. Avoid harsh or abrasive material that could scratch or in any way damage the optical surface.

7.2 STERILIZATION

The LED Rover is provided non-sterile. The device should be cleaned and disinfected or sterilized prior to use, according to your institution's requirements.

Caution *Do not reprocess the LED Rover using steam sterilization or autoclave. Use of these processes will result in damage to the instrument, and void its warranty.*

Process per your institution's validated procedures and process parameters.

Compatible Sterilization Method: STERRAD™ in STERRAD sterilizer. Double wrap with sterilization wrap. The full cycle of the Sterrad sterilizer comprises two identical exposure phases (6 minute injection, 2 minute diffusion and 2 minute plasma) preceded by a 10-minute pre-exposure plasma phase.

Follow the instructions provided with the STERRAD Sterilization System for sterilization processing via this equipment and method.

NOTE: The LED Rover has a black anodized exterior finish which is affected by Sterrad sterilization. The black anodized finish gradually bleaches; however, this is irrelevant to the functionality of the unit.

7.3 HIGH LEVEL DISINFECTION (IMMERSION)

The LED Rover may be immersed in high level disinfectant solution no deeper than six (6) inches (15cm) for as long as one (1) hour. Process according to the disinfectant manufacturer's instructions or your institution's established protocol. We make no claims other than the Rover's immersibility as stated.

8. LIMITED WARRANTY

Your fiberoptic products(s) carry a 1-year warranty from the date of shipment on workmanship and all defects of material except batteries. Should your product prove to have such defects within one year from the date of shipment, **ST TECHNOLOGIES** will repair or replace the product or component part without charge. Should your fiberoptic product(s) need servicing under this warranty, please contact **ST TECHNOLOGIES** for return authorization documentation. You should carefully pack product in a sturdy carton and ship it to the factory. Please include a note describing the defects, your name, telephone number and a return address. Warranty does not cover equipment subject to misuse, accidental damage, and normal wear and tear. This warranty gives you specific legal rights and you may also have other rights that vary from state to state.

- In the U.S call: 877 814 2237 (toll free)
- Customer Service 904 208 2291

9. CHART OF MEDICAL DEVICE SYMBOLS USED

	Manufacturer
	Date of manufacture (YYYYMM or YYYY)
	"Authorized Representative" in the European Community
	Caution, consult accompanying documents
	Consult Instructions for Use
	CE mark

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<p>CUDA[®] <i>SURGICAL</i></p> <p>ROVER III LED LIGHTSOURCE, ACMI PORT Caution: Federal law restricts this device to sale by or on the order of a licensed practitioner.</p> <p>REF LLS-3LBA LOT WOXXXXX</p> <p> ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.</p> <p> 20YY</p> <p>EC REP RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom</p> <p>R_x</p> <p>CE</p>	<p>CUDA[®] <i>SURGICAL</i></p> <p>ROVER III LED LIGHTSOURCE, ACMI PORT RECHARGEABLE BATTERIES Caution: Federal law restricts this device to sale by or on the order of a licensed practitioner.</p> <p>REF LLS-3LBA-R LOT WOXXXXX</p> <p> ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.</p> <p> 20YY</p> <p>EC REP RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom</p> <p>R_x</p> <p>CE</p>
<p>CUDA[®] <i>SURGICAL</i></p> <p>ROVER III LED LIGHTSOURCE, WOLF PORT Caution: Federal law restricts this device to sale by or on the order of a licensed practitioner.</p> <p>REF LLS-3LBW LOT WOXXXXX</p> <p> ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.</p> <p> 20YY</p> <p>EC REP RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom</p> <p>R_x</p> <p>CE</p>	<p>CUDA[®] <i>SURGICAL</i></p> <p>ROVER III LED LIGHTSOURCE, WOLF PORT RECHARGEABLE BATTERIES Caution: Federal law restricts this device to sale by or on the order of a licensed practitioner.</p> <p>REF LLS-3LBW-R LOT WOXXXXX</p> <p> ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.</p> <p> 20YY</p> <p>EC REP RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom</p> <p>R_x</p> <p>CE</p>

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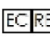
<p>CUDA[®] <i>SURGICAL</i></p> <p>ROVER III LED LIGHTSOURCE, ACMI PORT Caution: Federal law restricts this device to sale by or on the order of a licensed practitioner.</p> <p>REF LLS-3LBA LOT WOXXXXX</p> <p> ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.</p> <p> 20YY</p> <p> RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom</p> <p>R_x</p> <p>CE</p>	<p>CUDA[®] <i>SURGICAL</i></p> <p>ROVER III LED LIGHTSOURCE, ACMI PORT RECHARGEABLE BATTERIES Caution: Federal law restricts this device to sale by or on the order of a licensed practitioner.</p> <p>REF LLS-3LBA-R LOT WOXXXXX</p> <p> ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.</p> <p> 20YY</p> <p> RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom</p> <p>R_x</p> <p>CE</p>
<p>CUDA[®] <i>SURGICAL</i></p> <p>ROVER III LED LIGHTSOURCE, WOLF PORT Caution: Federal law restricts this device to sale by or on the order of a licensed practitioner.</p> <p>REF LLS-3LBW LOT WOXXXXX</p> <p> ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.</p> <p> 20YY</p> <p> RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom</p> <p>R_x</p> <p>CE</p>	<p>CUDA[®] <i>SURGICAL</i></p> <p>ROVER III LED LIGHTSOURCE, WOLF PORT RECHARGEABLE BATTERIES Caution: Federal law restricts this device to sale by or on the order of a licensed practitioner.</p> <p>REF LLS-3LBW-R LOT WOXXXXX</p> <p> ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.</p> <p> 20YY</p> <p> RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom</p> <p>R_x</p> <p>CE</p>

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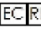
<p>CUDA[®] <i>SURGICAL</i></p> <p>ROVER III LED LIGHTSOURCE, ACMI PORT</p> <p>Caution: Federal law restricts this device to sale by or on the order of a licensed practitioner.</p> <p>REF LLS-3LBA LOT WOXXXXX R_x</p> <p> ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.</p> <p> 20YY CE</p> <p> RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom</p>	<p>CUDA[®] <i>SURGICAL</i></p> <p>ROVER III LED LIGHTSOURCE, ACMI PORT RECHARGEABLE BATTERIES</p> <p>Caution: Federal law restricts this device to sale by or on the order of a licensed practitioner.</p> <p>REF LLS-3LBA-R R_x LOT WOXXXXX</p> <p> ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.</p> <p> 20YY CE</p> <p> RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom</p>
<p>CUDA[®] <i>SURGICAL</i></p> <p>ROVER III LED LIGHTSOURCE, WOLF PORT</p> <p>Caution: Federal law restricts this device to sale by or on the order of a licensed practitioner.</p> <p>REF LLS-3LBW LOT WOXXXXX R_x</p> <p> ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.</p> <p> 20YY CE</p> <p> RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom</p>	<p>CUDA[®] <i>SURGICAL</i></p> <p>ROVER III LED LIGHTSOURCE, WOLF PORT RECHARGEABLE BATTERIES</p> <p>Caution: Federal law restricts this device to sale by or on the order of a licensed practitioner.</p> <p>REF LLS-3LBW-R R_x LOT WOXXXXX</p> <p> ST Technologies 6018 Bowdendale Ave. Jacksonville, FL 32216 U.S.A.</p> <p> 20YY CE</p> <p> RMS UK LTD 28 Trinity Rd. Nailsea, Somerset BS48 4NU United Kingdom</p>

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